



For Research Use


1. Connect the HDMI and USB sync cables to the laptop.

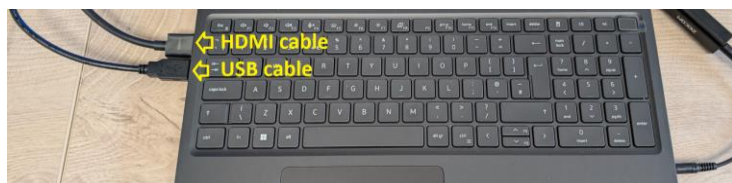


2. To duplicate the laptop screen to the display in the scanner room, press 'Windows' and 'P' keys, and select 'Duplicate'.

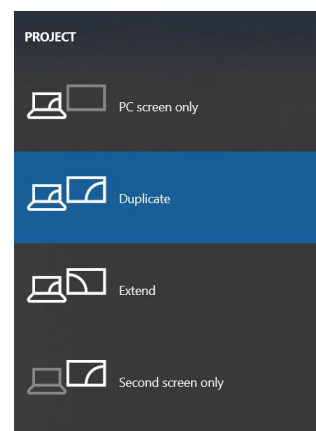


3. (Optional) If audio is required, take the audio cable from the UHS computer and connect to the laptop.

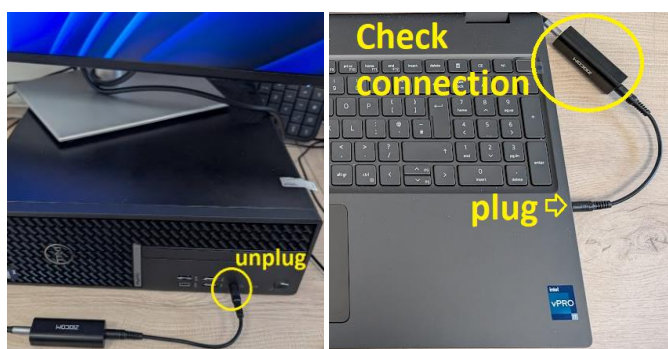

4. Now, check the SyncBox. By default, this should be displaying on the screen.

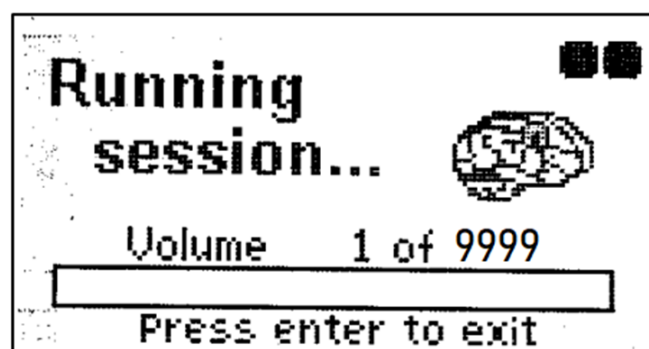

5. To test the triggers, scanner triggers will appear as letter 's'. ResponseGrips buttons correspond to the letters 'a', 'b', 'c', and 'd' for left thumb, left index, right index, and right thumb.

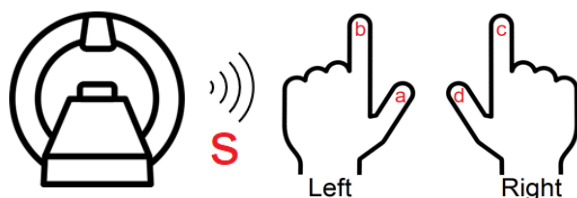

6. If all tests pass, load your paradigm now and make sure it's ready for the upcoming fMRI acquisitions.



Housekeeping: After all is finished, unplug the HDMI and USB, and re-connect the sound jack. Please leave everything else as is. **DO NOT** switch off the SyncBox and Fibre optic transmitter.

For queries, please contact MRI Physics (mriphysics@uhs.nhs.uk).